

# FINGER TIP PULSE Oximeter

## USER'S MANUAL

Model:  
 C101H1  C101A2  C101A3  
 C101B1  C101B2

Version: V1.2  
Release date: December 2020

**Responsibility of the Manufacturer**  
 IMDK only considers itself responsible for any effect on safety, reliability and performance of the equipment if: Assembly operations, extensions, re-adjustments, modifications or repairs are carried out by persons authorized by IMDK, and the electrical installation of the relevant room complies with national standards, and the instrument is used in accordance with the instructions for use. The equipment compliant with IEC60601-1 requirements of electrical safety and ensure the designated device's voltage and current meet the requirements of this Manual.

**Measuring principle**  
 Oximeter is based on the measuring principle haemoglobin, oxygenation of hemoglobin in the red and infrared light absorption characteristics in the region on the basis of the application "LambertBeer" Law of data presented. The instrument works by photoelectric detection of blood oxygen combined with the pulse volume recording technology, specific process is as follows:  
 First, the emission wavelength of used fluorescent tubes 660nm Red and wavelength 940nm Near-infrared light irradiation on the nails by photosensor measured signal. This information about wavelength range can be especially useful to clinicians.  
 To obtain data by electronic circuits and microprocessors, are displayed in Easy to read on. Operation schematic diagram:

**Section1 Safety**  
**1. Safety information**  
 Carefully read this manual about all safety informations, operation and specifications before using this oximeter.  
 •Do not place the equipment in children, pets and other places can be touched.  
 •This device is not intended for treatment, it can't be used to assess the accuracy of a pulse oximeter probe or a pulse oximeter monitor.  
 •Do not attempt to service the pulse oximeter. Only qualified service personnel should attempt any needed internal servicing.  
 •Do not use this device in situations where alarms required.  
 Although this oximeter provides alarm function, but the alarm does not meet IEC60601-1-8.  
 •The environment temperature should be guaranteed (working temperature: 5°C~40°C, transport and storage temperature: -10°C~+40°C). **When the ambient temperature is low or high, ensure that the product is recovered to room temperature before use.**  
 •It is not suitable for long-time continuous patient monitoring. Continual measurement must not exceed 2 hours. Do not charge during measurement. **Transfer of blood oxygen saturation and pulse rate data value oximetry in 8-10 seconds and data update cycle, more than 20 seconds.**  
 •SpO2 measurements may be adversely affected in the presence of high ambient light. Shield the sensor area (with a surgical towel, or direct sunlight, for example) if necessary.  
 •The following reasons will cause interference:  
 ■High-frequency electrocurgical  
 ■Placement of a sensor on an extremity with a blood pressure cuff arterial catheter, or intravascular line  
 ■The patient has hypotension severe vasoconstriction severe anemia or hypothermia  
 The patient is in cardiac arrest or is in shoe  
 ■Fingernail polish or false fingernails may cause inaccurate SpO2 readings

**2.1 Warnings**  
**WARNING: EXPLOSION HAZARD** Do not use the oximeter in a flammable atmosphere where concentrations of flammable anesthetics or other materials may occur.  
**WARNING:** The effects of degraded sensors and electrodes, or loosened electrodes, that can degrade performance or cause other problems.  
**WARNING:** Do not attempt to recharge normal dry-cell batteries, they may leak. And may cause a fire or even explode.  
**WARNING:** Do not use the pulse oximeter in an MRI or CT environment.  
**CAUTION:** Keep the operating environment free of dust, vibrations, corrosive, of flammable materials and extremes of temperature and humidity.  
**CAUTION:** The battery must be taken out from the battery compartment if the device will not be used for a long time.  
**CAUTION:** Do not operate the unit if it is damp or wet because of condensation or spills. Avoid using the equipment immediately after moving it from a cold environment to a warm, humid location.

**1.4 Definitions and Symbols**  
**1.4.1 Cleaning instruments:**  
 Use 70% of the alcohol cotton, clean silica gel sleeve, and test finger, probe and cavity cone. Please ensure that the instrument is inverted during cleaning to prevent the liquid from entering the instrument.  
 Don't put any liquid inside the instrument.  
 The instrument does not require the maintenance and calibration of the schedule except for the replacement of the battery.  
 Clinical test is a method commonly used to determine the oxygen accuracy. The measured arterial hemoglobin had an oxygen saturation, and the measurements were compared with the determined results of the arterial blood samples analyzed by the co-oximeter.  
 Measure ten times a day for ten minutes. It could take five years.  
 When plugging your finger into the Oximeter, your nail surface must be upward.

There are several reasons for a weak signal:  
 Low perfusion  
 The oximeter improper positioning  
 Cold temperatures and general health can cause low blood pressure  
 7. The pictures and interfaces in this manual are for reference only  
 8. This Manual is prepared based on the most complete configuration. Some configurations and functions may be not available in your

Declaration : Please use the medical alcohol to clean the rubber before each test and clean the tested finger with alcohol before and after the test. (The rubber inside of the Oximeter adopts medical rubber, which has no toxin, no harm, and brings no side effect such as allergy to the our skin).  
 Install two AAA batteries into battery cassette before covering its cover.  
 •Plug one finger into rubber hole of the Oximeter (it is best to plug the finger thoroughly) before releasing the clamp with the nail upwards.  
 •Press button on the front panel;  
 •Don't tremble your finger when the Oximeter is working. Your body is not recommended on moving status.  
 •Press the button on the front panel, if we want change display direction;  
 •Read relevant datum from display screen.  
 •If there is no signal input, oximeter can shut off automatically.  
 •Please replace new batteries when indicates the batteries are in low power.  
**Section 2 Introduction**  
**2.1 Brief Device Description**  
 This device intended for use in measuring and displaying functional oxygen saturation of arterial hemoglobin (%SpO2) and pulse rate (PR). The device measures SpO2 and PR with a SpO2 sensor and displays on the after certain further processing it can be used to measure human Hemoglobin Saturation and heart rate through finger  
**2.2 Application:**  
 The Oximeter intended to be used for home care and medical outpatient measure pulse oxygen saturation and pulse rate of adults. Oximetry is suitable for vascular diseases, respiratory diseases, the elderly over the age of 60, working more than 12 hours a day, extreme sports and alpine hypoxia environment of blood oxygen monitoring population and chronic alcoholism


**2.3 Symbols**

Symbols	Definition of symbols
	Follow instructions for us
	No alarm
	Type BF Applied Part
	Battery indication
<b>SpO2</b>	oxygen saturation of arterial hemoglobin
<b>PR</b>	Pulse Rate
	The symbol indicates that the device should be sent to the special agencies according to local regulations for separate collection after its useful life.
	The symbol indicates that the device complies with the European Council Directive 93/42/EEC concerning medical devices
	Authorized representative in the European community
	Serial Number
	Manufacturer
	Refer to use manual before application

**2.4 Specifications**

1. Type of protection against electric shock: Internally powered equipment  
 2. Degree of protection against electric shock: Type BF  
 3. Protection Against Ingress of Liquids: IP22 (protected against ingress of water when the water is dripping vertically and the monitor is tilted up to 15°)  
 4. Mode of operation: Continuous  
 5. Expected Service Life: 5 years  
 6. Display Type: more color display  
 7. SpO2:  
 Measurement range: 70%-99%  
 Accuracy: ±2% on the stage of 70%-99%, Unspecified (≤70%)  
 Resolution: ±1%  
 8. PR:  
 Measurement range: 30BPM-240BPM  
 Accuracy: ±2BPM or ±2%  
 9. Working Power:  
 Power Supply :2AAA 1.5V batteries ,  
 electric current: <50MA  
 Battery life: 2 AAA 1.5 V alkaline battery can be used for 30 hours continuously.  
 Battery voltage: Low battery indicator appears before battery power is lowered to normal operation  
 10. Dimension: 60\*35\*35mm  
 11. Environmental:

Operation Temperature	5°C-40°C
Storage Temperature	-10°C-40°C
Operation Humidity	15%-80%
Storage Humidity	15%-80%
Air Pressure	70-106kpa

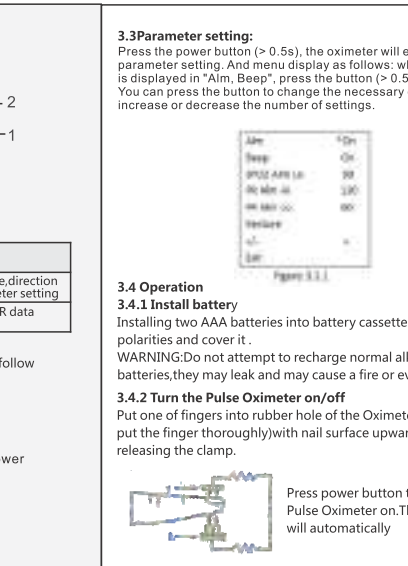

**3.3 Parameter setting:**  
 Press the power button (> 0.5s), the oximeter will enter into parameter setting. And menu display as follows: when the "" signal is displayed in "Alm, Beep", press the button (> 0.5s) to set on or off. You can press the button to change the necessary data. Select or increase or decrease the number of settings.

Figure 3.1.1

**3.4 Operation**  
**3.4.1 Install battery**  
 Installing two AAA batteries into battery cassette in correct polarities and cover it.  
**WARNING:** Do not attempt to recharge normal alkaline batteries, they may leak and may cause a fire or even explode.  
**3.4.2 Turn the Pulse Oximeter on/off**  
 Put one of fingers into rubber hole of the Oximeter (it is best to put the finger thoroughly) with nail surface upward, then releasing the clamp.

Press power button to turn the Pulse Oximeter on. The oximeter will automatically

be powered off when no finger in the device for longer than 16 seconds.  
**3.4.3 Read correspondent data from display screen.**  
**3.4.4 Display Description of**  
 The display interface of can rotate four directions with six different display modes after pressing the power button for less than 0.5s. It is shown as below:

**display modes**

Features display  
 •Two more color display, more display modes  
 •Low-power consumption, continuously four direction adjustable  
 •Low voltage indicator  
 •Automatically power off in 8 seconds when there is no signal  
 •Small in volume, light in weight and convenient to carry

**Section 4 maintenance and solution**  
**4.1 Maintenance and Preservation**  
 •Replace the batteries timely when low voltage lamp is on.  
 •Clean the surface of fingertip oximeter before it is used to diagnose patients.  
 •Remove the batteries inside if you will not operate the Oximeter for a long time.  
 •It would be better to preserve the product in -10~40 C (14-104 F) and humidity is 10%-80%.  
 •It is recommended that the product should be kept dry anytime. A wet ambience might affect its life time and even damage the product.  
 •Please follow the law of the local government to deal with used batteries.

**4.2 Product declaration**  
 Guidance and manufacture's declaration-electromagnetic radiation for other EQUIPMENTS and SYSTEMS

Radiation Test	Compliance	Electromagnetic environment-guidance
RF interference CISPR 11	Group 1	RF signal of Pulse Oximeter is simply created by its internal function. Therefore, its RF interference is very low and is not likely to cause any interference to nearby Electronic equipment
RF interference CISPR 11	Class B	The Pulse Oximeter applies to all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes

**4.3 Possible problems and effective solutions**

Problem	Possible reason	Solution
SpO2 or PR can not be shown normally	1. Finger is not plugged correctly 2. Patient's Oxymoglobin value is too low to be measured	1. Retry by plugging the finger 2. Try more times. If you can make sure there is no problem in the product, please go to hospital timely for exact diagnosis
SpO2 or PR is shown unsteady	1. The finger might not be plugged deep enough 2. Finger is trembling or the patient is on movement status	1. Retry by plugging the finger 2. Please remain at rest
The Oximeter can not be turned on	1. Inadequate power or power off 2. Batteries might be installed incorrectly 3. The Oximeter might be damaged	1. Please replace the batteries 2. Please reinstall the batteries 3. Please contact with local customer service center
Indication lamps are suddenly off	1. The product automatically shuts off when no signal is detected in 8 seconds. 2. Inadequate power	1. Normal 2. Replace the batteries

**4.4 Oximeter probe fault**  
**The fault of the probe is one of its own faults, and the other is caused by the external factors of the probe.**

1. If there is no SpO2 value in the oximeter, no red light emission will be detected by the probe. Poor probe. Replace the backup probe. SpO2 value returned to normal, determine the probe fault  
 2. The oximeter has no SpO2 value and has red light emission. It may be that the photo cell is insensitive to the light, the photoelectric tube is aging or the wire is broken Failure

3. probe external factors mainly in noise, jitter or patient finger keratinization, serious, probe launch, receiving part and unclean leakage. As long as the interference source is found, the probe is cleaned and the connecting mouth is properly used, the spot can be eliminated.

**Section 5. Applicable models**  
 C101H1, C101A2, C101A3, C101B1, C101B2.  
 Note: 1. The picture in the manual may be slightly different from the actual instrument.  
 2. Technical parameters are subject to change without prior notice.

**Section 6. Contact Information**  
 If you have any question about maintenance, technical specifications or malfunctions of devices, contact your local distributor or manufacturer.

-----  
 Shenzhen IMDK Medical Technology Co., Ltd.  
 904, 9F, Guangming Tianan Cloud Park Building, 255 Zhenmei Road, Zhenmei Community, Xihu Street, Guangming District, Shenzhen, PEOPLE'S REPUBLIC OF CHINA.  
 Post: 518107 Tel: +86-755-36637905

REC REP MedNet EC-REP GmbH, Borkstrasse 10, 48163, Münster, Germany